INSTRUCTIONS FOR USE

OneSett™

 Subcutaneous Administration Set

EN	OneSett™ Subcutaneous Administration Set – Instructions for use (page 2)
DA	OneSett™ Sæt til subkutan administration – Brugsanvisning (side 7)
FI	OneSett™ Subkutaaninen annostelusarja – Kayttoohjeet (sivu 12)
FR	OneSett™ Kit d'administration sous-cutanée – Model d'emploi (page 17)
NO	OneSett™ Subkutan administreringssett – Bruksanvisning (side 22)
SV	OneSett™ Subkutan administreringsset – Bruksanvisning (sidan 27)
TR	OneSett™ Subkutan Uygulama Seti – Kullanım talimatı (sayfa 32)
DE	OneSett™ subkutan Infusionsset – Gebrauchsanweisung (Seite 38)

	Caution	STERILE R	Sterilized using irradiation	Rx Only	Prescription only
	Consult instructions for use		Do not resterilize	MD	Medical device
	Do not use if package is damaged		Do not reuse		Sterile barrier system
	MR Unsafe		Temperature limit		Humidity limitation

INNOVATIVE HEALTH SCIENCES LLC

 1108 Kings Highway, Suite #4
 Chester, NY 10918 USA
 +1-855-680-0630

www.innohealthsci.com

04-1038 Rev 05 2026-08-14

EN INSTRUCTIONS FOR USE *Quick Start Guide***Indications for Using the OneSett™**

The OneSett™ is intended for use with the Insignis™ Syringe Infusion System in the subcutaneous administration of 13-17cP 20% immunoglobulin solutions with infusion systems operating with a safe 13.5psi constant pressure. Consult the drug manufacturer's package insert for recommendations.

Precautions and Warnings

- The administration sets are single use only. Reusing a single-use device could cause the device to not perform as intended or may result in infection. Do not re-sterilize.
- IHS advises against off-label use with the Insignis™ Syringe Infusion System.
- Never prime the infusion set or attempt to free a clogged tubing line while the set is inserted.
- Do not stretch or kink the delivery sets.
- It is recommended to store the administration sets in a dry place (humidity levels 15–90%) at room temperature (20–25°C; 68–77°F).
- Failure to store in correct storage conditions may lead to possible leak at the controller.

**Caution**

- Use of the OneSett™ is intended only for the patient for whom the device is prescribed.
- If using the infusion set for the first time, do so in the presence of a healthcare provider.
- Use recommended Becton Dickinson and Company® 50ml #300865 or similar syringes with the OneSett™.
- Before use, carefully inspect the administration set package. Do not use if the package has been opened or is damaged.
- Always follow the infusion instructions located on the package insert of the manufactured drug.
- The Insignis™ Syringe Infusion System has not been tested in a Magnetic Resonance Imaging (MRI) environment and therefore cannot be recommended for use in such an environment.
- Verify the drug volume in the syringe(s) prior to beginning the infusion.
- Store infusion sets in a cool, dry place.
- Do not leave infusion sets in direct sunlight or inside a vehicle.
- If for any reason you experience a medical emergency during use, call your local emergency number and/or your healthcare provider immediately.
- Use of flow rate calculation software and/or computation documentation (tables, charts, graphs, etc.) may be unnecessary for use with the OneSett™ as the flow dial indicates the exact flow rate delivered to the patient.
- Always consult and comply with guidance provided by your healthcare providers and drug manufacturers regarding the location and the number of infusion sites.
- Each OneSett is calibrated from 10-60ml/hr/site for the number of needles included. Each needle site can deliver up to 60ml/hr of 13-17cP 20% immunoglobulin solution. A four-site needle configuration can provide a total flow rate of 240ml/hr (60ml/hr/site).

INSTRUCTIONS FOR USE

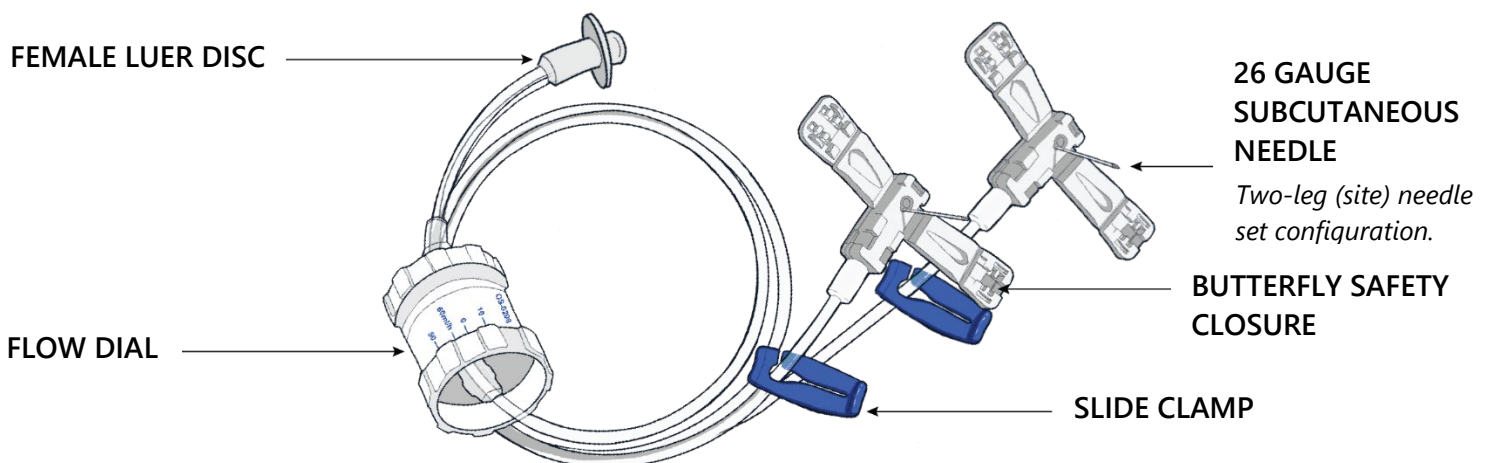
EN **Caution**

- If a needle site is clamped-off while using a multi-needle set, the flow rate to the remaining needle site(s) will increase. To compensate for this increase, it is advised to manually decrease the flow rate on the flow dial. See flow rate compensation recommendations page 6.
- When using a two, three, or four needle set configuration, over-infusion may occur to remaining sites if one or more of the needles is blocked.
- Inaccurate medication delivery, infection, and/or site irritation may result from improper needle insertion and maintenance of the infusion site.
- Ensure that the needle is at 90°.
- Be sure to remove the needle guard prior to insertion.
- Since individuals inherently have varying amounts of subcutaneous fat on the body, be sure to select the appropriate needle length to fit your needs.
- For suspected site reactions, perform subsequent infusions in alternative sites.
- To stop the flow immediately, use the slide clamp and/or set the OneSett™ dial to 0ml/hr.
- Ensure the syringe is empty before removing the syringe from the driver. If medication is still in the syringe, activate the lever once to infuse the remaining medication.

Adverse Event Reporting

Any serious incident that has occurred in relation to the device should be reported to the manufacturer and the competent authority of the Member State in which the user and/or patient is established. Contact information for both is shown on the last page of this document.

OneSett™ Subcutaneous Administration Set (*Selectable Rate Flow Dial and Needle Administration Set – 24" tubing*):

EN **INSTRUCTIONS FOR USE** *Quick Start Guide***Starting the Infusion for Subcutaneous Administration using the OneSett™:****1. Ensure you have all supplies for your infusion.**

Verify that you are using the correct administration set, making sure to check the number of sites and needle length on the packaging.

**2. Wash hands**

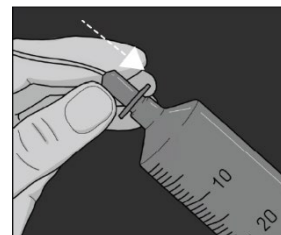
Wash hands thoroughly and if needed, wear disposable gloves.

**3. Fill Syringe**

Make sure the drug product is at room temperature and fill one or more 50ml BD® or similar syringe(s) to your required dose. Note: it may take several syringes to receive your prescribed dose.

**4. Attach OneSett™**

Using aseptic technique, remove the end cap on OneSett™ and attach to the luer on 50ml BD® or similar syringe; twist to lock.

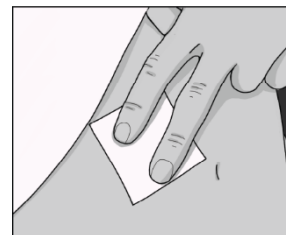
**5. Prime OneSett™**

Set OneSett™ to 60ml/hr. Gently press on syringe plunger and prime set, taking care to stop before medication reaches the needles. Set OneSett™ to 0ml/hr (OFF).

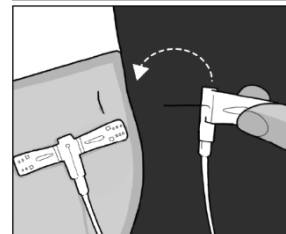
EN **Starting the Infusion for Subcutaneous Administration** *using the OneSett™:***6. Prep Infusion Site(s)**

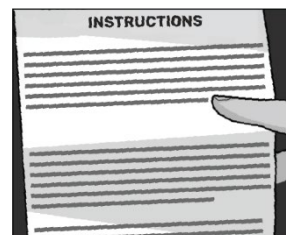
Using aseptic technique, prep infusion site as directed by clinician. Remove rubber band from butterfly, fold wings back, and remove needle guard.

**7. Insert Needle(s)**

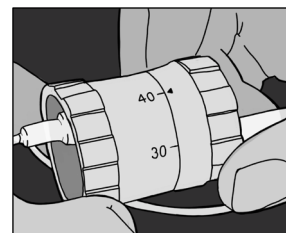
Insert needle(s) into appropriate locations in subcutaneous tissue at 90° angle. Lay butterfly wings flat on skin. Cover site with tape or adhesive dressing.

**8. Activate Infusion Pump**

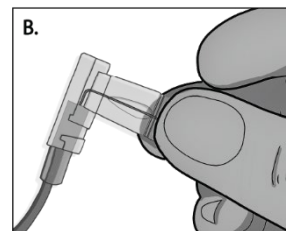
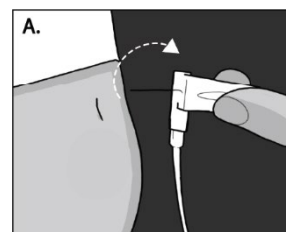
Refer to the infusion pump instructions for use provided by the manufacturer to load and activate the infusion pump.

**9. Set Flow Rate**

As advised by your healthcare provider, set appropriate flow rate and begin administration. Do not exceed the maximum flow rate as indicated by your provider.

**Ending the Infusion****1. Remove Dressing**

Pause before removing needles; allow site pressures to decrease naturally. Remove needle dressing. To minimize accidental risk of needle stick injury, use one hand to remove needle. Using the same hand, squeeze butterfly safety closure wings together to firmly encase needle.

**2. Dispose of Materials**

Carefully dispose of materials in a sharps container per local regulations.

INSTRUCTIONS FOR USE

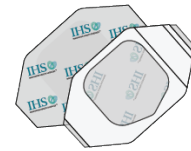
EN Product Identifiers

OneSett™ Subcutaneous Administration Set (box of 10)

Configuration	6mm	9mm	12mm	14mm
1 site	OS-0106	OS-0109	OS-0112	OS-0114
2 sites	OS-0206	OS-0209	OS-0212	OS-0214
3 sites	OS-0306	OS-0309	OS-0312	OS-0314
4 sites	OS-0406	OS-0409	OS-0412	OS-0414

Sterile Dressing

Sterile dressing/adhesive is included in the administration set box.

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Definition of Symbols

	Caution	REF	Catalog number	SN	Serial number
	Consult instructions for use	LOT	Batch code	MD	Medical device
	Use by YYYY-MM-DD	QTY	Quantity	UDI	Unique device identifier
	Date of manufacture	STERILE R	Sterilized using irradiation	Rx Only	Prescription only
	Temperature limit		Do not resterilize		Do not reuse
	MR Unsafe		Humidity limitation		Do not use if package is damaged
	Sterile barrier system	EU REP	European Authorized Representative		Manufacturer
CE	European Conformity	CH REP	Swiss Authorized Representative		Importer

CH REP

OncoMedical AG
 Hintere Hauptgasse 9
 CH-4800 Zofingen,
 Switzerland

EU REP

MedEnvoy Global B.V.
 Prinses Margrietplantsoen 33 –
 Suite 123 2595 AM The Hague
 The Netherlands

UKRP

MedEnvoy UK Limited
 167–169 Great Portland Street
 5th Floor
 London, W1W 5PF
 United Kingdom

INNOVATIVE HEALTH SCIENCES LLC

1108 Kings Highway, Suite #4
 Chester, NY 10918 USA
 +1-855-680-0630
 www.innohealthsci.com

CE 2797

OneSett™ Subcutaneous Administration Set is a trademark of Innovative Health Sciences LLC (IHS).
 IHS is ISO 13485 certified and abides by the quality policies as outlined in 21 CFR 820.